



disposable. versatile. cost-effective.

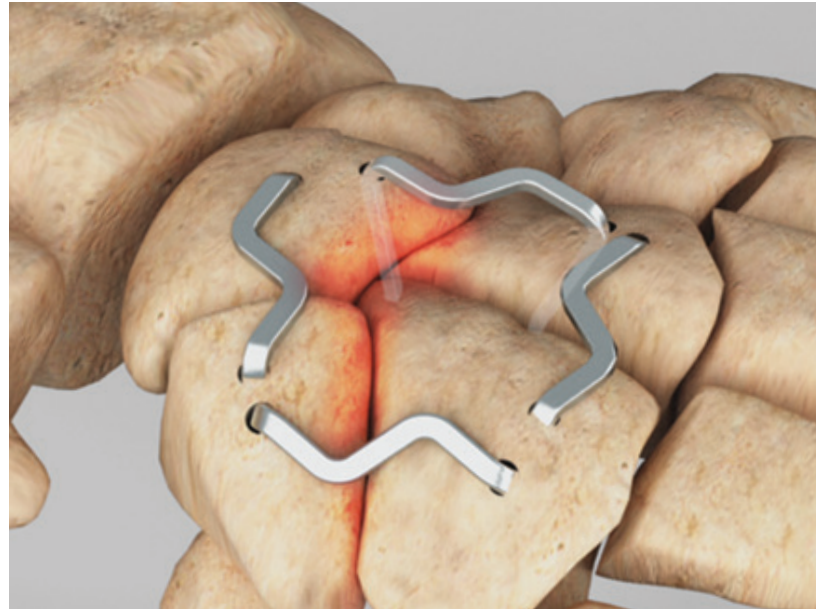

A SIMPLE STAPLE FIXATION SYSTEM...

for arthrodesis, osteotomies, and other skeletal fixation procedures. The InstaFix shape memory fixation system is body temperature activated. The use of Nitinol, a biocompatible shape memory alloy, allows dynamic compression of the InstaFix staple implant upon reaching body heat. A S-Back design provides continuous and active bi-cortical compression at the surgical site.

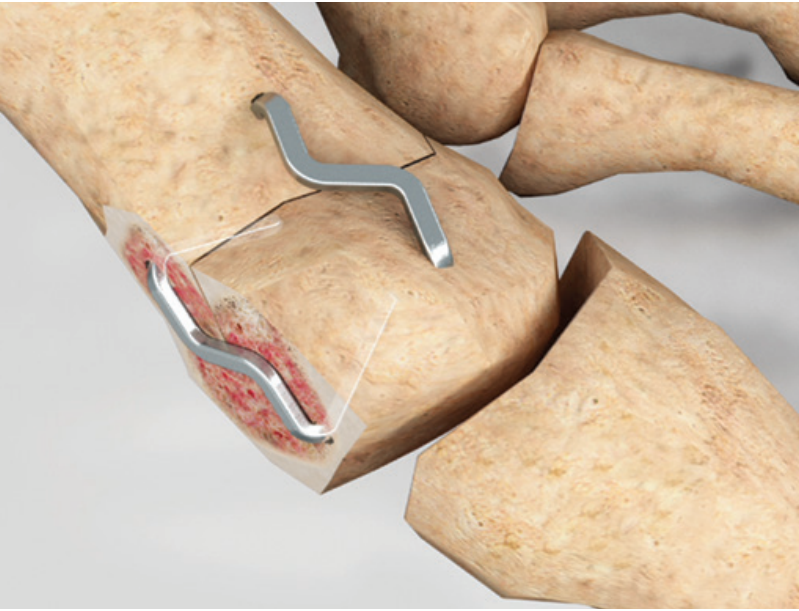

Implant & Accessory Kit

Sterile & fully disposable

- Cost savings for hospital or surgery center
- Eliminates pre & post-operative set management

Body temperature activated Nitinol implant

- Instantaneous and continuous compression
- Eliminates the need for external heating
- Simple alternative to plates & screws

Ergonomically designed instrumentation

- Reduces surgical complexity
- Provides maximum visualization of surgical site


Low profile implant design

- Helps minimize tissue disruption
- Offers unobstructed views for post-op X-rays




InstaFix™ in the open position

Bi-cortical compression when staple is implanted

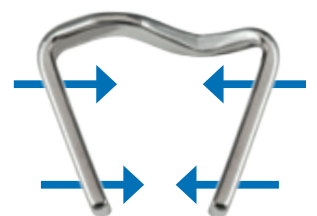

ORDERING INFORMATION

InstaFix™ Implant Kit

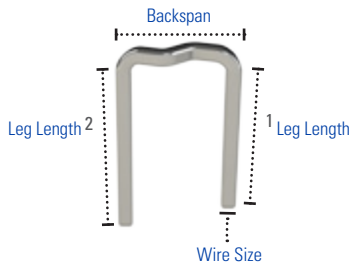

Accessory Kit

Catalog No.	Description	Kit Size	Use with InstaFix™ implant Kit
400-1115-SP	InstaFix™ Accessory Kit*	1.5 mm	400-3701-SP
400-1120-SP	InstaFix™ Accessory Kit*	2.0 mm	400-3703-SP, 400-3709-SP, 400-3715-SP, 400-3721-SP
400-1135-SP	InstaFix™ Accessory Kit*	3.5 mm	400-3722-SP, 400-3723-SP, 400-3724-SP, 400-3725-SP

Catalog No.	Description	Backspan	Leg Length	Wire Size	Drill Size	Use with Accessory Kit
400-3701-SP	InstaFix™ Implant Kit	8 mm	8 mm	1.2 mm	1.5 mm	400-1115-SP
400-3703-SP	InstaFix™ Implant Kit	10 mm	10 mm	1.5 mm	2.0 mm	400-1120-SP
400-3709-SP	InstaFix™ Implant Kit	15 mm	10 mm	1.5 mm	2.0 mm	400-1120-SP
400-3715-SP	InstaFix™ Implant Kit	10 mm	15/17 mm	1.5 mm	2.0 mm	400-1120-SP
400-3721-SP	InstaFix™ Implant Kit	15 mm	15/17 mm	1.5 mm	2.0 mm	400-1120-SP
400-3722-SP	InstaFix™ Implant Kit	20 mm	15 mm	3 mm x 2 mm	3.5 mm	400-1135-SP
400-3723-SP	InstaFix™ Implant Kit	20 mm	20 mm	3 mm x 2 mm	3.5 mm	400-1135-SP
400-3724-SP	InstaFix™ Implant Kit	25 mm	20 mm	3 mm x 2 mm	3.5 mm	400-1135-SP
400-3725-SP	InstaFix™ Implant Kit	30 mm	20 mm	3 mm x 2 mm	3.5 mm	400-1135-SP



Optional Size Kit

Catalog No.	Description
400-7777-SP	InstaFix™ Sizers Kit

*Accessory kit contains (1) drill & (2) locating pins

MULTIPLE SIZES FOR MULTIPLE INDICATIONS

The InstaFix shape memory fixation staple system is versatile, offering nine distinct sizes for all extremities procedures.

